

# Certificate of Analysis


## Retatrutide 20mg

**Compound** : Retatrutide  
**Lot number** : 2025-04-07  
**Analysis date** : 2025-05-01  
**Purity %** : 99.35%  
**Quantity** : 23.86 mg  
**Method** : HPLC-UV-MS

**Client** : Transforma Peptides  
[www.transformapeptides.com](http://www.transformapeptides.com)

CAS number : 2381089-83-2

### High Performance Liquid Chromatography (HPLC) UV – Purity and Quantity Test

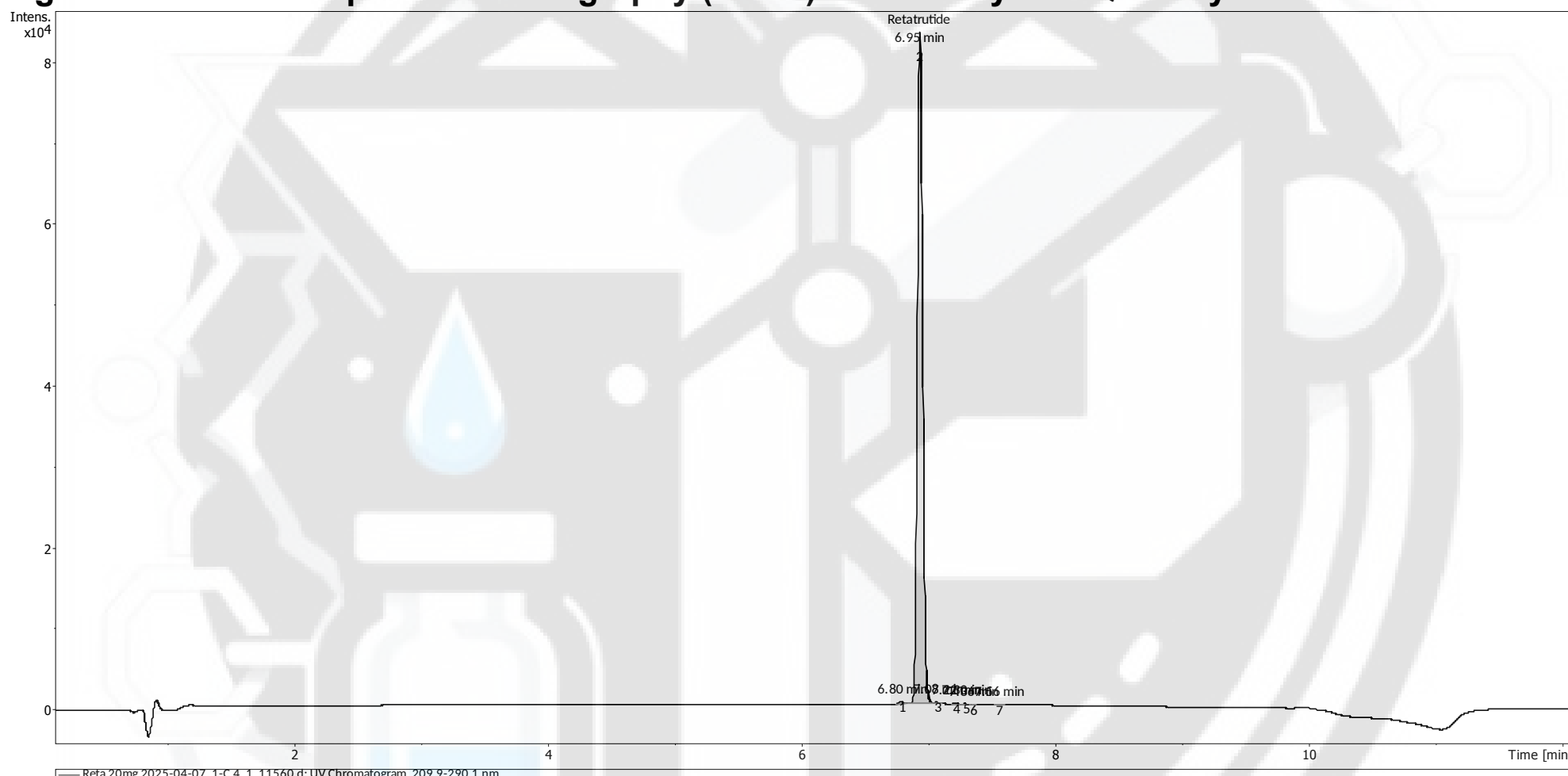

| PEAK LIST |            | Number of detected peaks: 7 |       |
|-----------|------------|-----------------------------|-------|
|           | Time (min) | Area                        | %Area |
| 1         | 6.80       | 7.59E+02                    | 0.31  |
| 2         | 6.95       | 2.44E+05                    | 99.35 |
| 3         | 7.08       | 4.42E+02                    | 0.18  |
| 4         | 7.22       | 2.30E+02                    | 0.09  |
| 5         | 7.30       | 7.16E+01                    | 0.03  |
| 6         | 7.36       | 4.87E+01                    | 0.02  |
| 7         | 7.56       | 5.15E+01                    | 0.02  |

Retatrutide

Analysis Performed by  
Ken Pendarvis, ChE  
Analytical Chemist  
MZ Biolabs  
[contact@mzbiolabs.com](mailto:contact@mzbiolabs.com)

**Quantification by HPLC-UV**  
Measured quantity : 23.86 mg/vial

Note: Injectable peptides may contain salts and sugars to aid in solubility and act as pH buffers.  
These are not normally detected using UV and are not considered impurities.



2025-05-07

# Retatrutide 20mg

CAS number : 2381089-83-2

## Mass Spectrometry (MS) – Identity Test

### Identity confirmed using HPLC-MS

Molecular weight calculated using monoisotopic m/z values from mass spectrum

Expected monoisotopic mass : 4728.47 Da


Measured monoisotopic mass : 4728.52 Da

**Molecular weight confirmed**

Note : Monoisotopic m/z values are not easily seen in full spectrum view for larger molecules and peptides.

The dominant isotopic peak (base peak) shown in the spectrum below can be used to approximate the average molecular weight frequently reported by vendors and databases as a secondary means of confirmation.

### Recorded MS spectrum

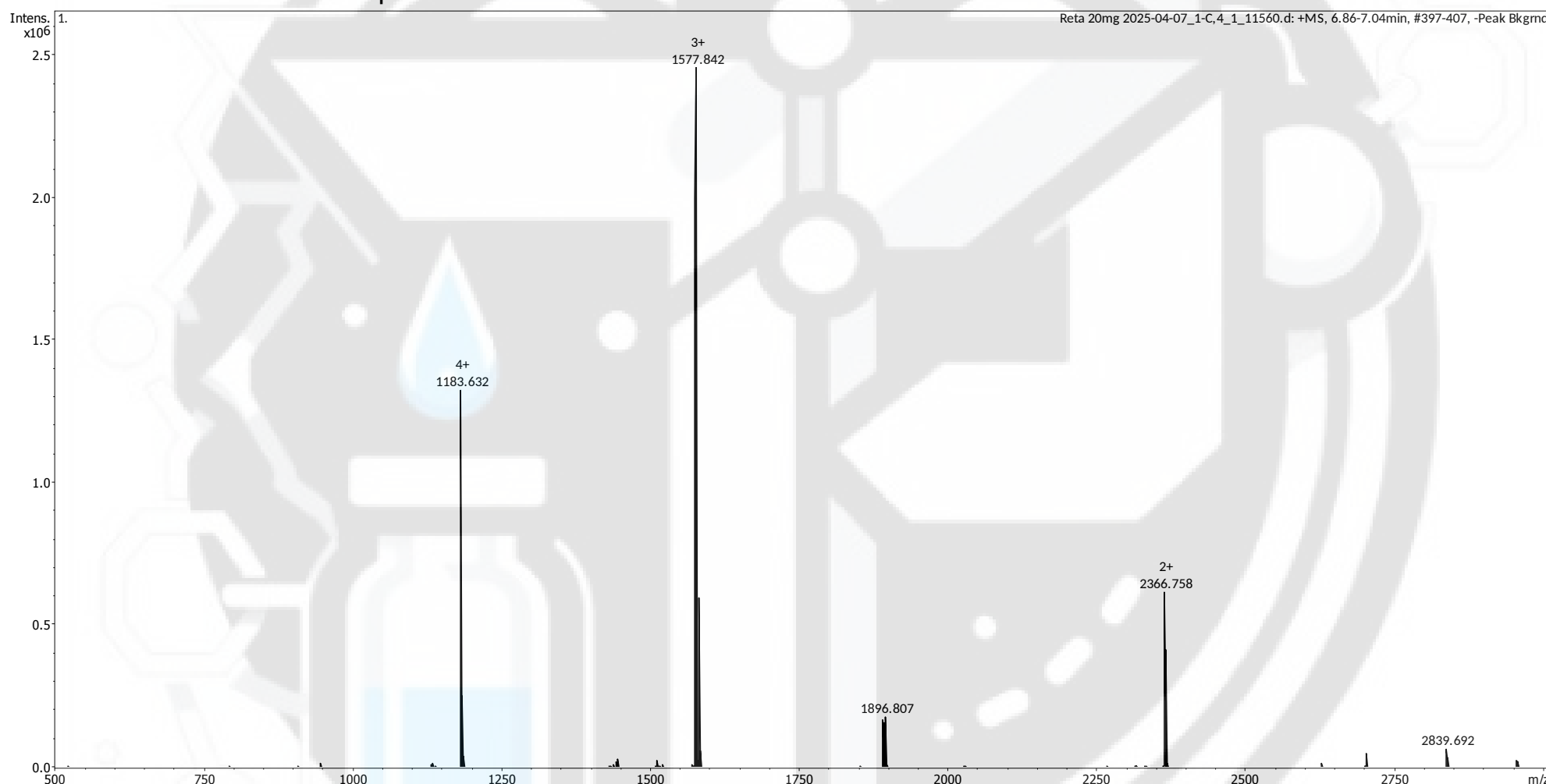

Analysis Performed by  
Ken Pendarvis, ChE  
Analytical Chemist  
MZ Biolabs  
[contact@mzbiolabs.com](mailto:contact@mzbiolabs.com)



2025-05-07